



Quantus™ Fluorometer Specifications

Size	22.7 × 11.5 × 4.5cm
Weight	400 grams
Linear Dynamic Range	Five (5) orders of magnitude, assay-dependent
Limit of Detection	1pg/ml double-stranded DNA with QuantiFluor® dsDNA Dye (50pg/ml using standard protocol) 3.9pg/ml RNA using QuantiFluor® RNA Dye
Detectors	Solid-state silicon sensor
Wavelengths	Excitation Filters: Red, 640nm shortpass; Blue, 495nm shortpass Emission Filters: Red, 660–720nm; Blue, 510–580nm
Calibration Type	Single-point (blank and standard)
Sample Tube	0.5ml PCR tubes (We recommend Axygen (Cat.# PCR-05-P)
Automatic Sleep Mode	After 1 hour of inactivity
Power Requirements	5V, 0.2A maximum
Computer Requirements (for optional Quantus™ Software)	Windows® 8, Windows® 7 (32 or 64 bit) .NET Framework 4.0 or later At least one USB port available
Warranty	1 year

Ordering Information

Quantus™ Fluorometer	Cat.# E6150
QuantiFluor® dsDNA System	Cat.# E2670
QuantiFluor® ssDNA System	Cat.# E3190
QuantiFluor® RNA System	Cat.# E3310

Products may be covered by pending or issued patents or may have certain limitations. Please visit our Web site for more information.
Helix, QuantiFluor and Quantus are trademarks of Promega Corporation
NanoDrop is a registered trademark of Thermo Fisher Scientific. PicoGreen is a registered trademark of Molecular Probes, Inc.
Qubit is a registered trademark of Life Technologies Corporation. Windows is a registered trademark of Microsoft Corporation.



PROMEGA CORPORATION • 2800 WOODS HOLLOW ROAD • MADISON, WI 53711-5399 USA • TELEPHONE 608-274-4330
www.promega.com • © 2013 PROMEGA CORPORATION • ALL RIGHTS RESERVED • PRICES AND SPECIFICATIONS SUBJECT TO CHANGE WITHOUT PRIOR NOTICE • PRINTED IN USA 06/13 • 22226-4807 • PART #BR249



Quantus
FLUOROMETER

An affordable, sensitive single-tube fluorometer for fluorescence-based quantitation of DNA and RNA.



Find Out How Much Nucleic Acid You *Really* Have.

Sensitivity

Use of the Quantus™ Fluorometer with QuantiFluor® Dyes provides sensitivity up to 40,000 times higher than that of absorbance-based quantitation using the NanoDrop® spectrophotometer, depending on the assay used.

Cost Savings

- Realize savings of up to 25% over similar instruments.
- No computer is required.
- No need for special accessories.

Better Results

Successful and accurate quantitation leads to higher success rates in downstream assays such as PCR, cloning, transfection and next-generation sequencing.

Convenience

Your new fluorometer might be as close as your nearest Helix™ on-site stocking cabinet.

Quantus™ Software

The Quantus™ Software is easy to install to allow data transfer from the Quantus™ Fluorometer to a computer. With the Quantus™ Software, you no longer have to record fluorescence measurements manually, and you can easily assign and change sample names using the computer keyboard.

Dual Channels

Measure fluorescence using two channels:

Blue fluorescence channel:

Excitation wavelengths up to 495nm (495nm shortpass), emission at 510–580nm

Red fluorescence channel:

Excitation wavelengths up to 640nm (640nm shortpass), emission at 660–720nm

Seamless Workflow

The Quantus™ Fluorometer is preprogrammed with protocols for Promega QuantiFluor® Dye Systems for single-stranded DNA, double-stranded DNA and RNA quantitation.

Flexibility

The Quantus™ Fluorometer is optimized for use with QuantiFluor® Dyes, but the fluorometer also is compatible with other fluorescent quantitation assays.

Specificity

Fluorescence-based nucleic acid quantitation methods are more specific than absorbance-based methods that rely on A₂₆₀ or A₂₈₀ measurement.

Ease of Use

Intuitive menus and a large LCD allow easy navigation.

The single-point calibration process saves time by eliminating the need to generate a standard curve for each set of unknown samples. The Quantus™ Fluorometer saves the calibration data.

Bench Space

The compact instrument is only 22.7 × 11.5 × 4.5cm.

10 TIMES MORE Sensitive Than Qubit® 2.0 Fluorometer

22.7cm

4.5cm

11.5cm

Simple & Flexible Workflow

Easy to Use

Optimized with Promega QuantiFluor® Dyes


Up to 25% SAVINGS Over Similar Instruments

Quantitate with Confidence.

High Sensitivity

Fluorescence-based nucleic acid quantitation methods are more sensitive than absorbance-based methods. This allows you to complement your absorbance-based methods and quantify your target molecule over a wider range of concentrations.

The high sensitivity makes the Quantus™ Fluorometer perfect when working with sample types such as formalin-fixed paraffin-embedded (FFPE) tissue that contain low levels of nucleic acid. You can confidently proceed to your downstream detection assays, knowing that there is DNA or RNA present.

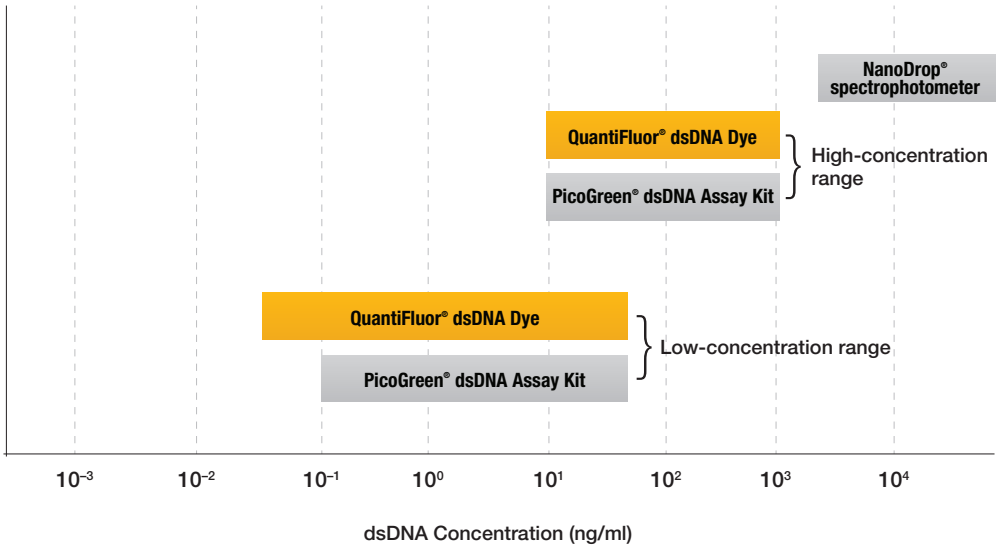

Superior Specificity

Fluorescence-based nucleic acid quantitation methods also are more specific than absorbance-based quantitation methods that rely on A₂₆₀ or A₂₈₀ measurement:

- Single-stranded DNA, double-stranded DNA and RNA absorb strongly at 260nm and cannot be distinguished by absorbance-based quantitation methods.
- Contaminants that absorb at 260nm, such as proteins or guanidine isothiocyanate, can cause an overestimation of nucleic acid concentrations.

